

Preliminary

Instruction Manual

MSX-Exxxx and LabVIEW 8.5



Product information

This manual contains the technical installation and important instructions for correct commissioning and usage, as well as production information according to the current status before printing.

The content of this manual and the technical product data may be changed without prior notice.

ADDI-DATA GmbH reserves the right to make changes to the technical data and the materials included herein.

Warranty and liability

The user is not permitted to make changes to the product beyond the intended use, or to interfere with the product in any other way.

ADDI-DATA shall not be liable for obvious printing and phrasing errors. In addition, ADDI DATA, if legally permissible, shall not be liable for personal injury or damage to materials caused by improper installation and/or commissioning of the board by the user or improper use, for example, if the board is operated despite faulty safety and protection devices, or if notes in the operating instructions regarding transport, storage, installation, commissioning, operation, thresholds, etc. are not taken into consideration. Liability is further excluded if the operator changes the board or the source code files without authorisation and/or if the operator is guilty of not monitoring the permanent operational capability of working parts and this has led to damage.

Copyright

This manual, which is intended for the operator and its staff only, is protected by copyright.

Duplication of the information contained in the operating instructions and of any other product information, or disclosure of this information for use by third parties, is not permitted, unless this right has been granted by the product licence issued. Non-compliance with this could lead to civil and criminal proceedings.

ADDI-DATA software product licence

Please read this licence carefully before using the standard software. The customer is only granted the right to use this software if he/she agrees with the conditions of this licence.

The software must only be used to set up the ADDI-DATA boards.

Reproduction of the software is forbidden (except for back-up and for exchange of faulty data carriers). Disassembly, decompilation, decryption and reverse engineering of the software are forbidden. This licence and the software may be transferred to a third party if this party has acquired a board by purchase, has agreed to all the conditions in this licence contract and the original owner does not keep any copies of the software.

Trademarks

- ADDI-DATA, MSX-Box und MSX-E sind eingetragene Warenzeichen der ADDI-DATA GmbH.
- Turbo Pascal, Delphi, Borland C, Borland C++ sind eingetragene Warenzeichen von Borland Software Corporation.
- Microsoft .NET, Microsoft C, Visual C++, Windows XP, Windows 98, Windows 2000, Windows 95, Windows NT, Windows Embedded NT, Windows Vista, Windows Server 2003, Windows Embedded, Windows Server 2000 und MS-DOS sind eingetragene Warenzeichen von Microsoft Corporation.
- LabVIEW, LabWindows/CVI, DasyLab, Diadem sind eingetragene Warenzeichen von National Instruments Corporation.
- CompactPCI ist ein eingetragenes Warenzeichen von PCI Industrial Computer Manufacturers Group.
- VxWorks ist ein eingetragenes Warenzeichen von Wind River Systems, Inc.
- RTX ist ein eingetragenes Warenzeichen von Ardence.

Warning

The following risks result from improper implementation and from use of the product contrary to the regulations:



Personal injury



Damage to the module, PC and peripherals



Pollution of the environment

■ **Protect yourself, the others and the environment!**

■ **Read carefully the safety precautions (yellow leaflet).**

If this leaflet is not with the documentation, please contact us and ask for it.

■ **Observe the instructions of the manual.**

Make sure that you do not forget or skip any step. We are not liable for damages resulting from a wrong use of the product.

■ **Used symbols**



IMPORTANT!

Designates hints and other useful information.



WARNING!

It designates a possibly dangerous situation.

If the instructions are ignored the module, PC and/or peripheral may be **destroyed**.



WARNING!

It designates a possibly dangerous situation.

If the instructions are ignored the module, PC and/or peripheral may be destroyed and persons may be injured or killed.

Contents

1	Introduction.....	5
2	How to import the Web Service of MSX-Exxxx in LabVIEW 8.5	6
3	How to use the Web Service in LabVIEW 8.5	14
4	How to receive data from the socket.....	17
5	Contact and support	18

Figures

Fig. 2-1:	„Diagram“ window	6
Fig. 2-2:	“Import Web Service”: WSDL URL	7
Fig. 2-3:	“Import Web Service”: User name and password	8
Fig. 2-4:	“Import Web Service”: Library and destination directory	9
Fig. 2-5:	“Import Web Service”: Methods to import.....	10
Fig. 2-6:	“Import Web Service”: Summary	11
Fig. 2-7:	“Import Web Service” Creating a Wrapper Library.....	12
Fig. 2-8:	“Import Web Service”: Completed	13
Fig. 3-1:	“Web Service”: User library	14
Fig. 3-2:	Open the Web Service	15
Fig. 3-3:	Using the other functions: Example	15
Fig. 3-4:	Closing the Web Service	16
Fig. 4-1:	Opening a TCP socket connection	17
Fig. 4-2:	Reading values from the socket.....	17
Fig. 4-3:	Closing the TCP connection.....	17

1 Introduction


This document describes the following:

- How to import the Web Service of an Ethernet module MSX-Exxxx in LabVIEW 8.5
- How to use the Web Service
- How to receive data from the socket

2 How to import the Web Service of MSX-Exxxx in LabVIEW 8.5


- Start LabVIEW and go to the diagram window.
- Go to "Tools\Import\Web Service".

Fig. 2-1: „Diagram” window

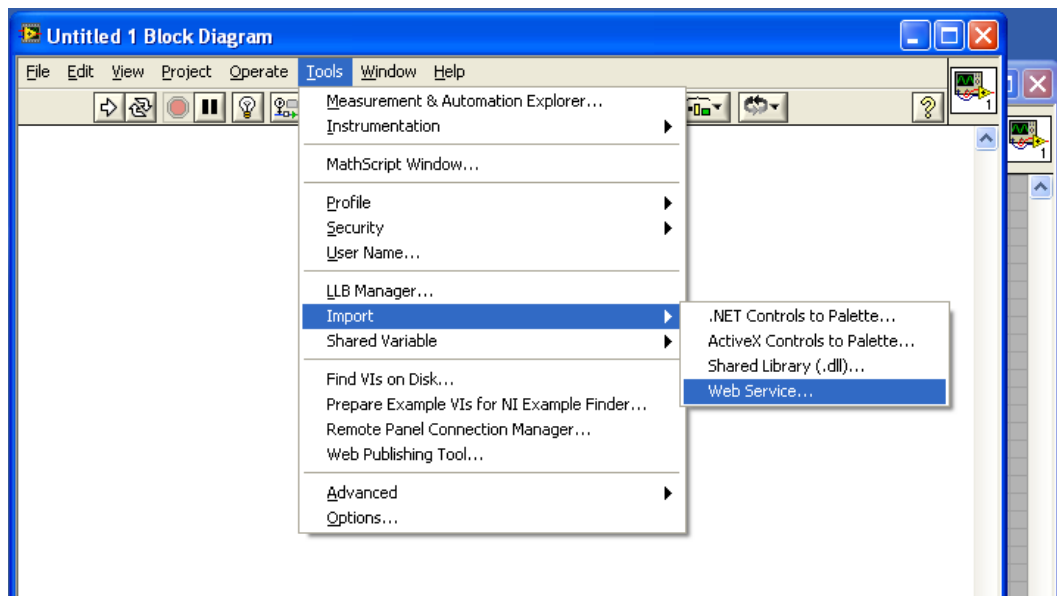

The first "Import Web Service" window appears.

Fig. 2-2: "Import Web Service": WSDL URL

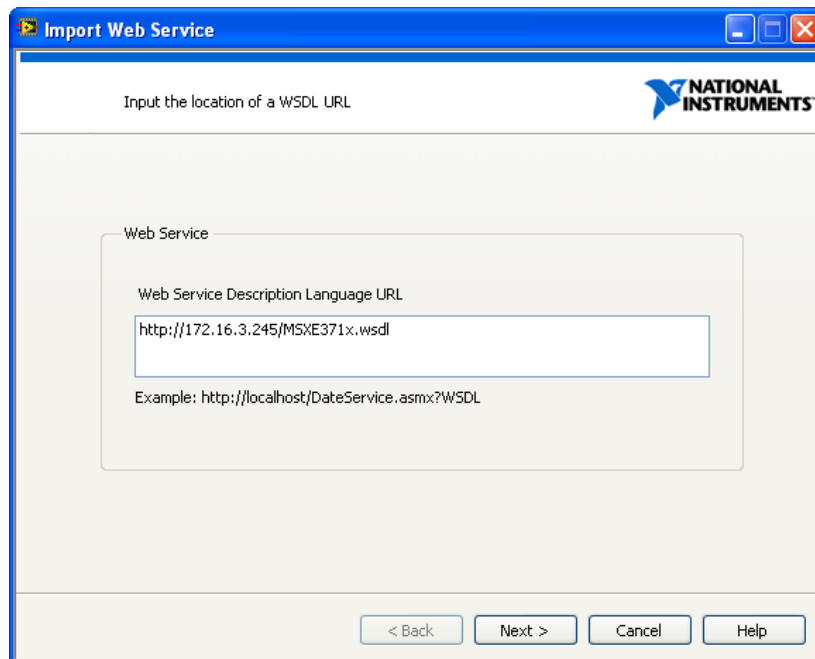

- Enter the Path of the WSDL file.
- Click on "Next".

The next "Import Web Service" window appears.

Fig. 2-3: “Import Web Service”: User name and password

Import Web Service

Input the authentication information

Authentication

☒ Web service requires authentication

User name* mxadmin Password *****

Domain

Proxy Server

☐ Use a proxy server for your LAN

Name/Address* Port 0

☐ Proxy service requires authentication


User name* Password

Domain

* Marked items are required when the option is enabled.


< Back Next > Cancel Help

- Enter the user name (mxadmin) and password (mxadmin) to access the web server.
- Click on “Next”.

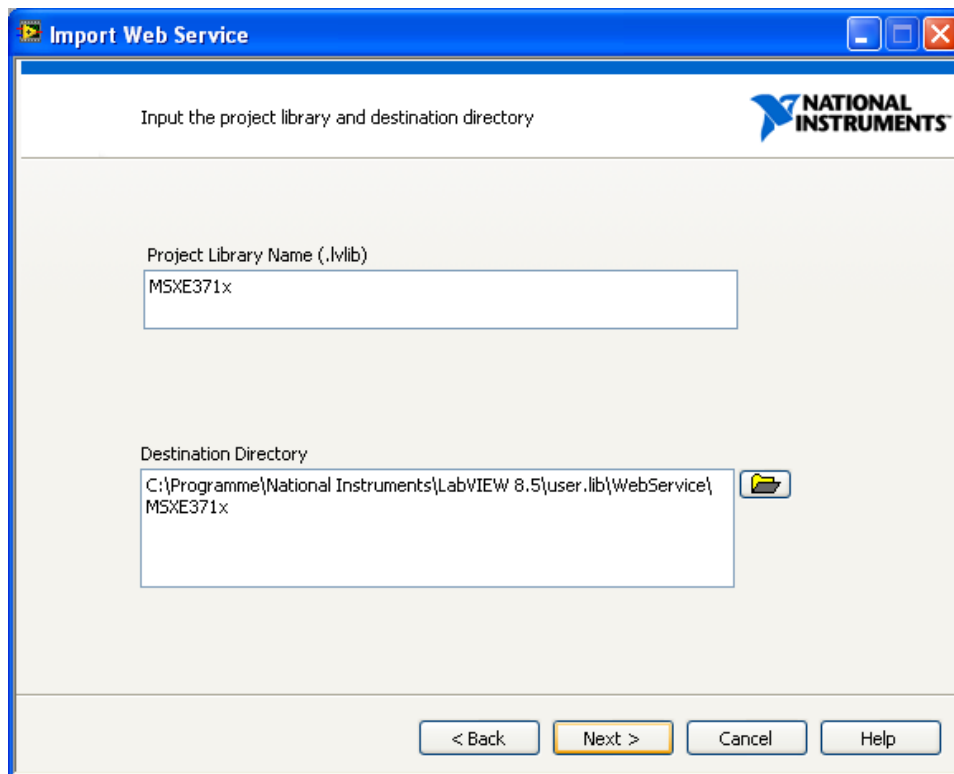Fig. 2-4: “Import Web Service”: Library and destination directory


- Click on “Next”.

A window to select the method you want to import appears.

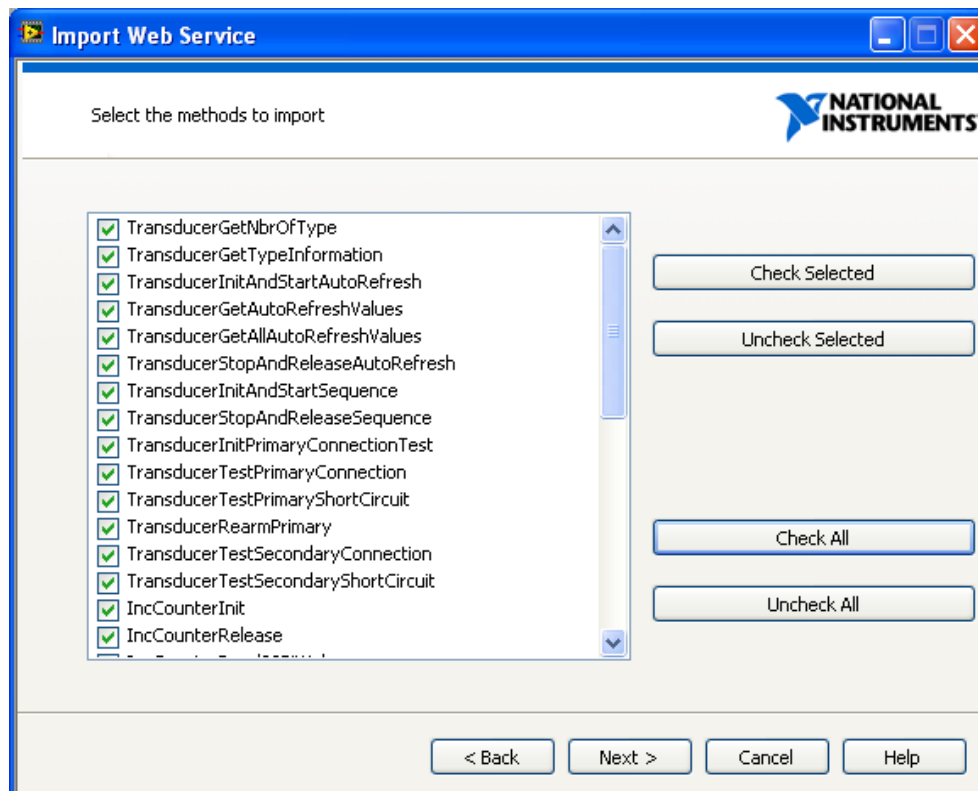Fig. 2-5: “Import Web Service”: Methods to import


- Select the method (=function) you want to import (we selected all the functions).
- Click on "Next".

Fig. 2-6: "Import Web Service": Summary

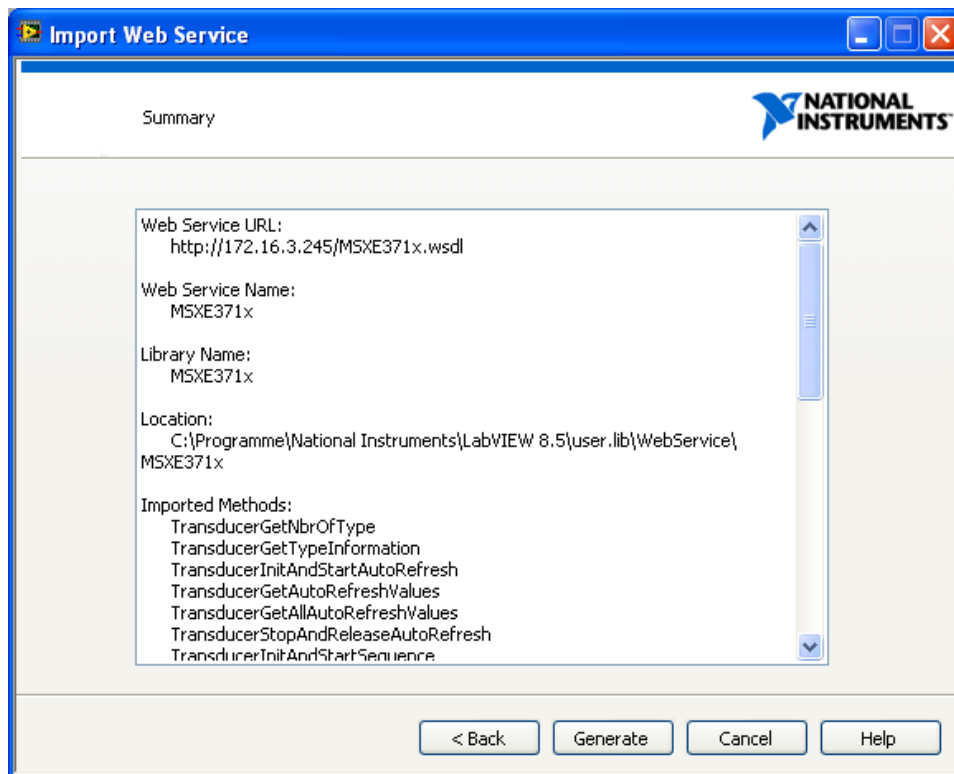


- Click on "Generate"

LabVIEW generates a Wrapper Library with the desired functions.

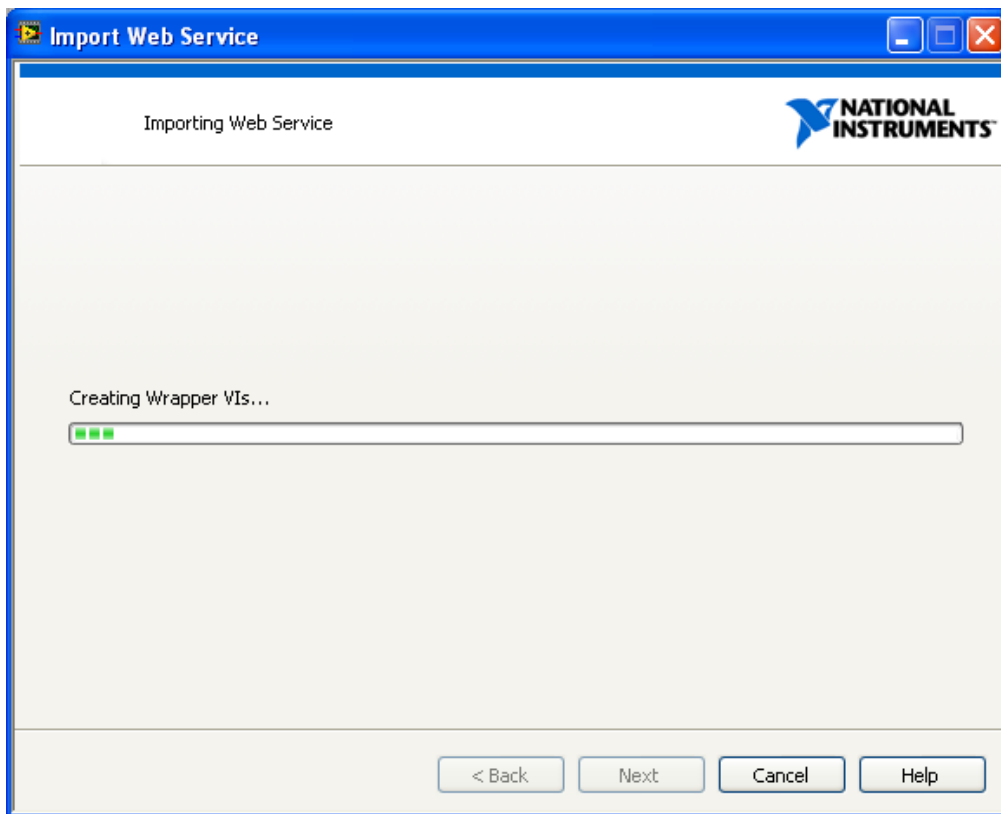Fig. 2-7: “Import Web Service” Creating a Wrapper Library

As soon as the import is completed, the following window appears.

Fig. 2-8: "Import Web Service": Completed

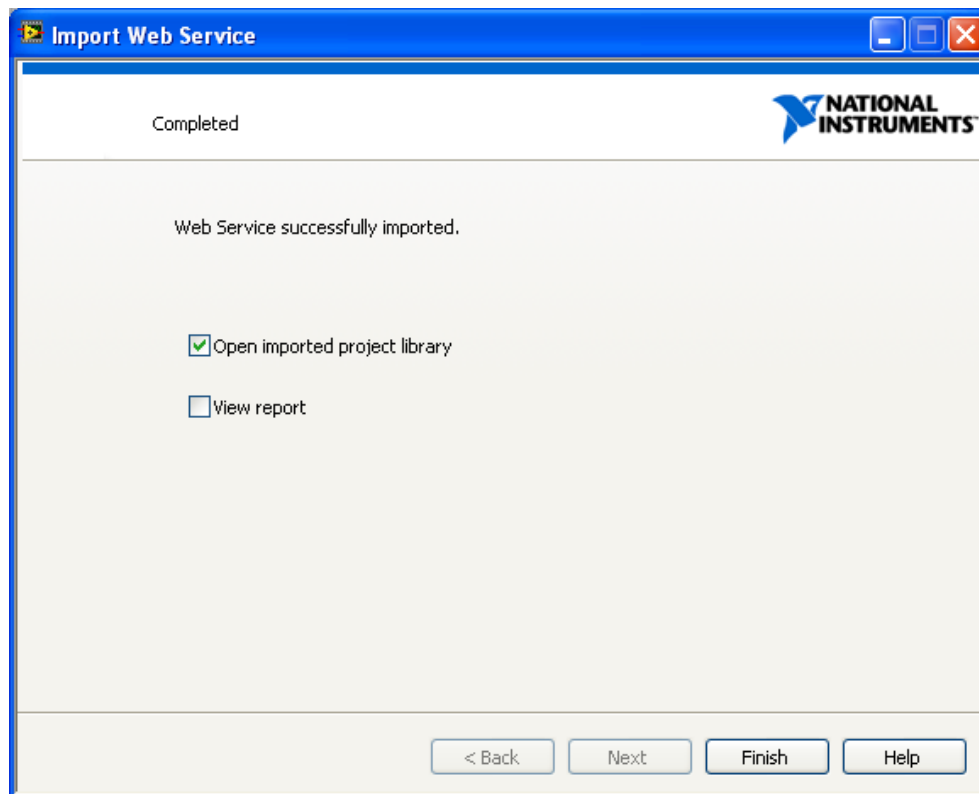


- Click on "Finish".
- Restart LabVIEW to be able to use the Web Service.

3 How to use the Web Service in LabVIEW 8.5


After a restart of LabVIEW you will find the new Web Service functions in the user library:

Fig. 3-1: “Web Service”: User library

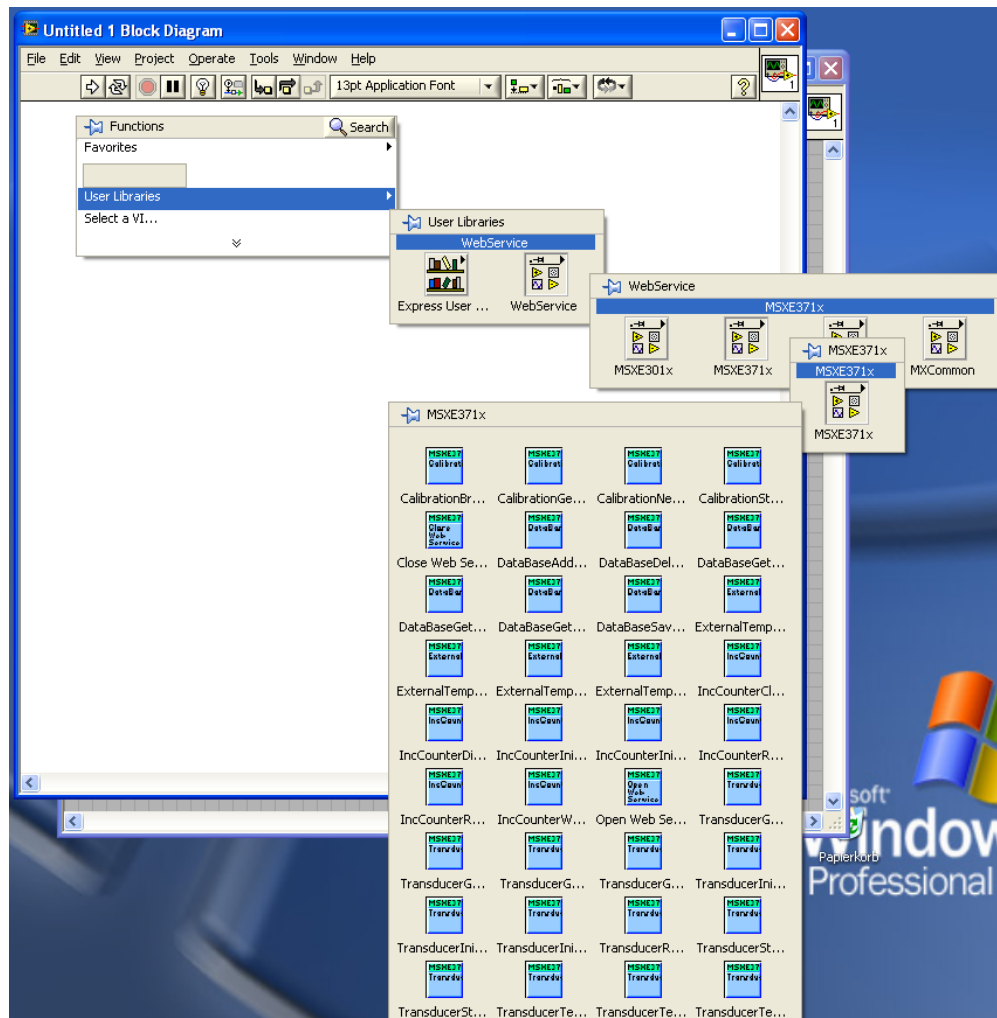


If you want to use it:

- At first, open the Web Service.

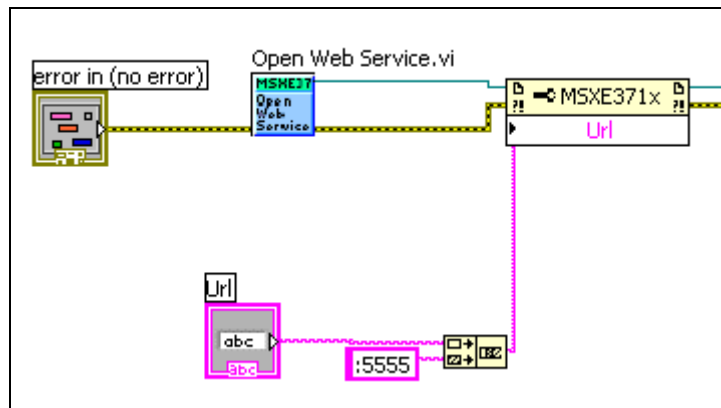
Fig. 3-2: Open the Web Service

Then you can use the other functions:

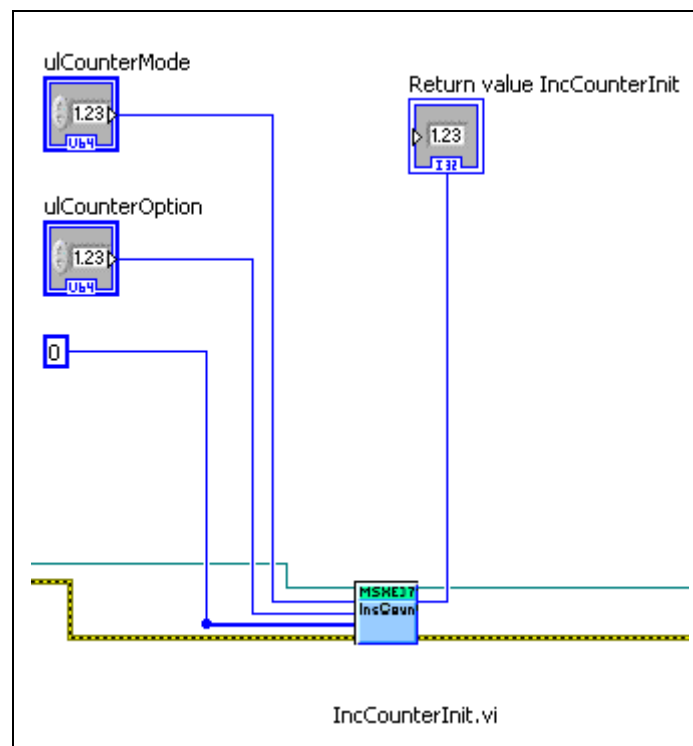Fig. 3-3: Using the other functions: Example

At the end of the application:


- Do not forget to close the Web Service

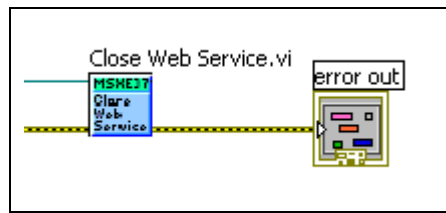Fig. 3-4: Closing the Web Service

4 How to receive data from the socket


- Open a TCP socket connection

Fig. 4-1: Opening a TCP socket connection

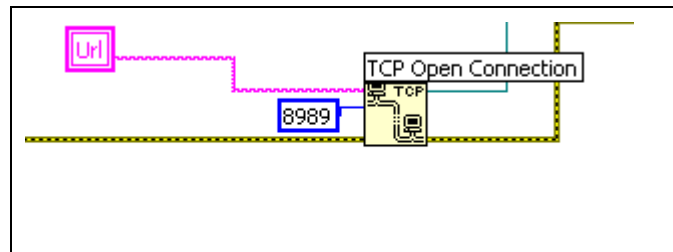


Now you can read the values from the socket:

Fig. 4-2: Reading values from the socket

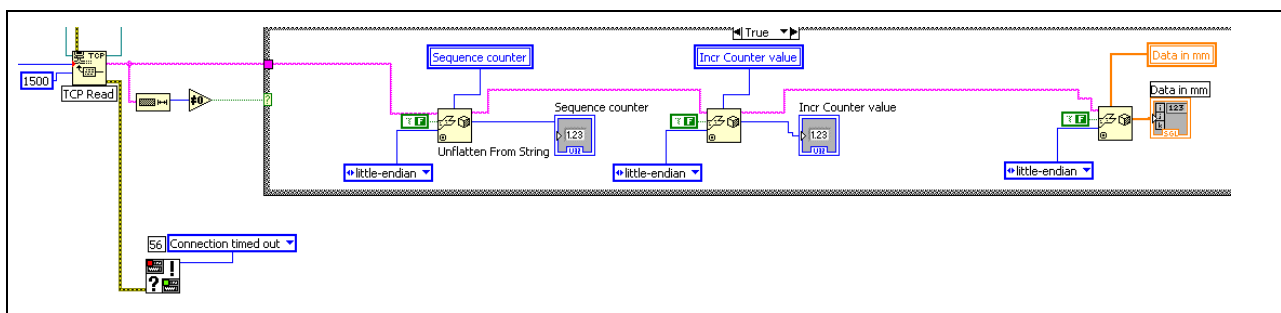

- Do not forget to close the TCP connection.

Fig. 4-3: Closing the TCP connection

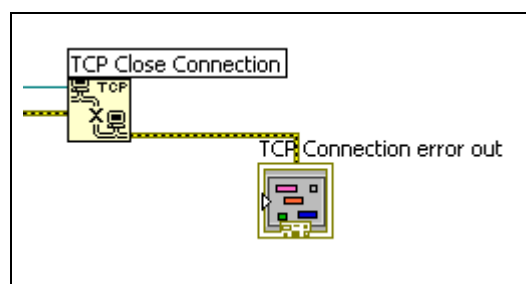

5 Contact and support

Do you have any questions? Write to us or phone us:

Address: ADDI-DATA GmbH
Airpark Business Center
Airport Boulevard B210
77836 Rheinmünster

E-Mail: info@addi-data.com

Phone: +49 7229 1847-0

Manual and software download on the internet:

www.addi-data.com